


Negative pressure table mounted
PU-1818 Power Unit captures particulates
and exhausts HEPA filtered air out of enclosure.

Ceiling Mounted PU-2525 Power Unit
delivers positive pressure
HEPA filtered air to Cybio Unit.

Positive pressure air
moves past Cybio instrument

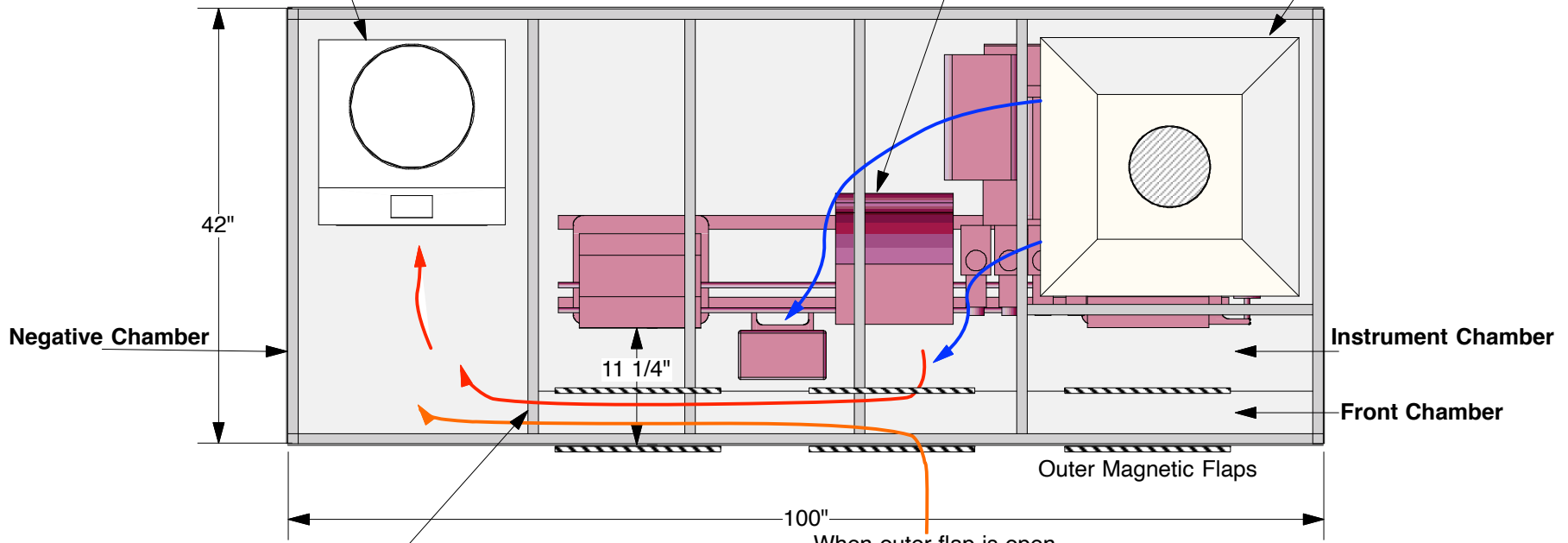

When outer flap is open,
ambient room air is drawn inward towards
the Negative Chamber without entering Instrument Chamber.

Elongated cutout in vinyl wall allows
positive pressure air to move into negative chamber.

Cybio Unit 67 1/2" x 32"

Cybio - SELMA 1536 Enclosure
Product & Personnel Protection
Enclosure OD: 100" W x 42" D x 52" H

Drawing Is Accepted By:
Signature:
Date:

BIO BUBBLE	
CONTROLLED ENVIRONMENTS CUSTOM SOLUTIONS	
Phone: 970.224.4262 Fax: 970.224.2419 Email: sales@biobubble.com	
Version #:	01
Date:	14 August 2015
Scale:	3/4" = 1' 0"
Drawn By:	Jack